


**EDIBON PATENT**




Always included in the supply:

Teaching Technique used

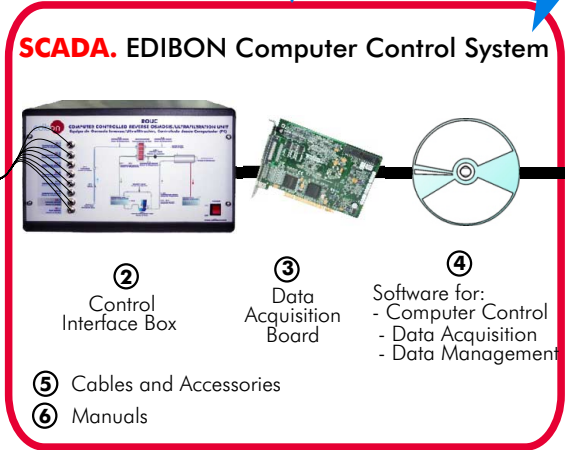

**OPEN CONTROL + MULTICONTROL + REAL TIME CONTROL**

① Unit: ROUC. Reverse Osmosis/Ultrafiltration Unit

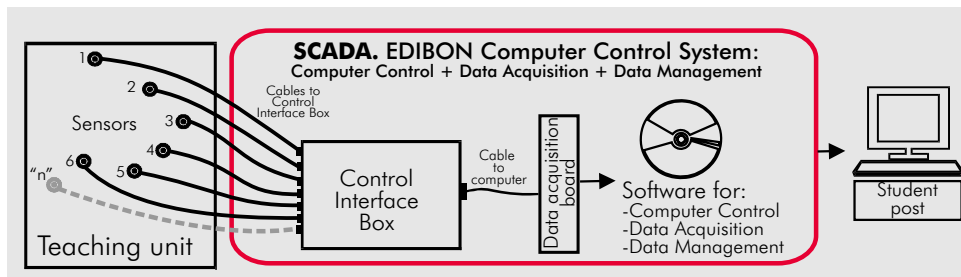

[www.edibon.com](http://www.edibon.com)  
 Products  
 Products range  
 Units  
 12.-Food & Water Technologies

## DESCRIPTION

Laboratory Scale Computer Controlled Reverse Osmosis/Ultrafiltration Unit designed to study and to provide practical training in these downstream processing techniques.

As well as the processing of the whey, the membranes can be also used to demonstrate, for example: clarification and concentration of fruit juices, potabilization of water, pre-concentration of milk for cheese manufacturing, etc.

This unit is composed by a three head positive displacement pump. The pump motor is controlled by an inverter that allows obtaining different flows depending on the required separation process. Because it is a positive displacement pump, the obtained flow is proportional to the motor speed.

Therefore, the pump supplies the effluent to the two tubular membranes connected in series placed in the membrane module. This solution keeps moving in the inner part of the tube, the permeate flows from the inner part of the tube towards the outer one and the concentrate is obtained in the final extreme of the tubes.

The ROUC unit makes a filtration process by cross flow, which avoids the appearance of impurities that reduces the capacity of the system.


The effluent moved by the pump receives some heat from the pump so its temperature increases. The permeate is stored in a final product tank. The concentrate goes back to the initial tank in order to be used again as an effluent but, due to the increase of temperature, this concentrate passes through a heat plate exchanger.

You have to know two different variables of the effluent in order to do the separation practices through Reverse Osmosis and Ultrafiltration:

- Chemical composition of the fluid.
- Characteristics of the membrane to know which one you should use in the process you want to obtain: the Retention Character, which is defined as the allowed pass size and the Flow type that, in this case, consists on cross flow but it depends on the membrane system that has been used and its measurement depends on the required production.

This Computer Controlled Unit, is supplied with the EDIBON Computer Control System (SCADA), including: Control Interface Box + Data Acquisition Board + Computer Control and Data Acquisition Software, for controlling the process and the parameters involved.

## PROCESS DIAGRAM AND ELEMENTS ALLOCATION

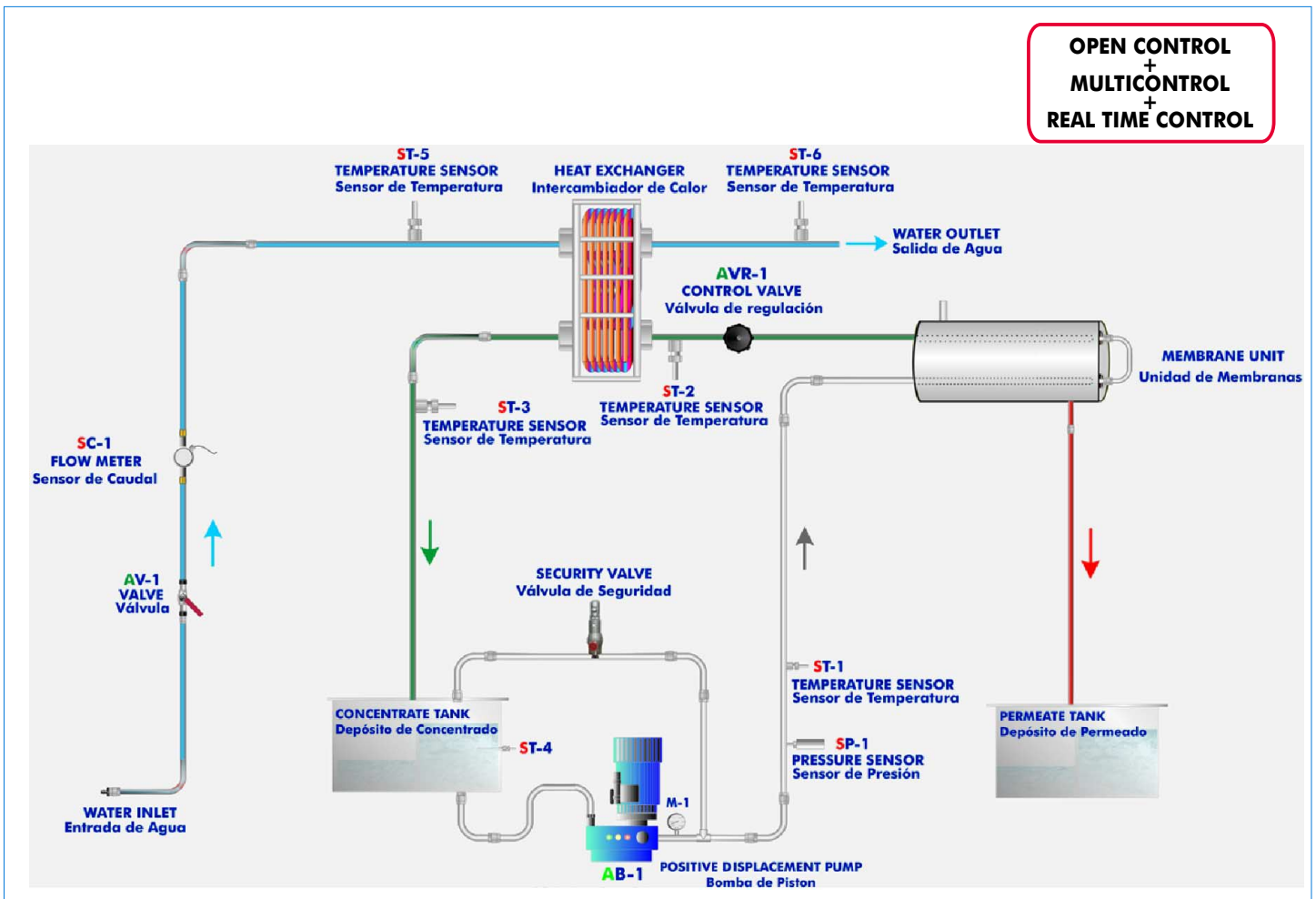

## SPECIFICATIONS

### Items supplied as standard

#### ① ROUC. Unit:

Anodized aluminium structure and panels in painted steel (epoxi paint).

Main metallic elements in stainless steel.

Diagram in the front panel with similar distribution to the elements in the real unit.

Stainless steel feed tank (15 l. approx. capacity).

Three head positive displacement feed pump:

Maximum flow: 38 l./min,

Maximum pressure: 150 bar, maximum operation pressure: 55 bar.

This pump is used to feed product to the membrane module.

Flow and pressure adjustable. Variable speed control.

The pump has a relief valve to protect the unit.

Inverter that controls the pump motor.

Membrane module:

Two tubular membranes connected in series.

Membrane diameter: 12.5 mm.

Membrane area: 0.024 m<sup>2</sup>.

Tubeside volume: 75 ml.

Maximum operation pressure: 55 bar.

Process control valve.

Plates heat exchanger for the concentrate.

Permeate stainless steel collecting tank (15 l. approx. capacity).

6 Temperature sensors, "J" type, range: -40 to 750°C.

Pressure sensor (0-100 bar).

Flow sensor (water inlet), range: 0-6 l/min.

2 Membranes of each model of the following are supplied:

RO01: Retention Character (99% NaCl), Process: Reverse Osmosis, Max. Pressure (45 bar).

UF02: Retention Character (20,000 Da), Process: Ultrafiltration, Max. Pressure (10 bar).

UF03: Retention Character (200,000 Da), Process: Ultrafiltration, Max. Pressure (15 bar).

The filtration membranes differ only in length and quantity from those installed in a real plant.

Rapid changeover from Reverse Osmosis to Ultrafiltration and back.

Easy cleaned and maintenance.

All elements in contact with the process fluid are constructed from hygienic design materials such as stainless steel, PTFE and silicone rubber.

The unit incorporates wheels for its mobility.

#### ② ROUC/CIB. Control Interface Box :

Control interface box with process diagram in the front panel and with the same distribution that the different elements located in the unit, for an easy understanding by the student.

All sensors, with their respective signals, are properly manipulated for -10V. to +10V computer output. Sensors connectors in the interface have different pines numbers (from 2 to 16), to avoid connection errors. Single cable between the control interface box and computer.

The unit control elements are permanently computer controlled, without necessity of changes or connections during the whole process test procedure.

Simultaneously visualization in the computer of all parameters involved in the process. Calibration of all sensors involved in the process. Real time curves representation about system responses. Storage of all the process data and results in a file. Graphic representation, in real time, of all the process/system responses.

All the actuators' values can be changed at any time from the keyboard allowing the analysis about curves and responses of the whole process. All the actuators and sensors values and their responses are placed in only one computer screen.

Shield and filtered signals to avoid external interferences.

Real time computer control with flexibility of modifications from the computer keyboard of the parameters, at any moment during the process. Real time computer control for pumps, compressors, resistances, control valves, etc. Open control allowing modifications, at any time and in a real time, of parameters involved in the process simultaneously.

Three safety levels, one mechanical in the unit, other electronic in control interface and the third one in the control software.

#### ③ DAB. Data Acquisition Board:

PCI Data acquisition board (National Instruments) to be placed in a computer slot. Bus PCI.

Analog input: Number of channels= 16 single-ended or 8 differential. Resolution= 16 bits, 1 in 65536.

Sampling rate up to: 250 KS/s (Kilo samples per second). Input range (V)= ±10V.

Data transfers=DMA, interrupts, programmed I/O. Number of DMA channels=6.

Analog output: Number of channels=2. Resolution= 16 bits, 1 in 65536. Maximum output rate up to: 833 KS/s.

Output range(V)= ±10V.Data transfers=DMA, interrupts, programmed I/O.

Digital Input/Output: Number of channels=24 inputs/outputs. DO or DI Sample Clock frequency: 0 to 1 MHz.

Timing: Counter/timers=2. Resolution: Counter/timers: 32 bits.

#### ④ ROUC/CCSOF. Computer Control+Data Acquisition+Data Management Software:

Compatible with actual Windows operating systems. Graphic and intuitive simulation of the process in screen.

Compatible with the industry standards. Registration and visualization of all process variables in an automatic and simultaneously way. Flexible, open and multicontrol software, developed with actual windows graphic systems, acting simultaneously on all process parameters.

Management, processing, comparison and storage of data. Sampling velocity up to 250,000 data per second guaranteed. Calibration system for the sensors involved in the process.

It allows the registration of the alarms state and the graphic representation in real time. Comparative analysis of the obtained data, after the process and modification of the conditions during the process. Open software, allowing to the teacher to modify texts, instructions. Teacher's and student's passwords to facilitate the teacher's control on the student, and allowing the access at different work levels.

This unit allows that the 30 students of the classroom can visualize simultaneously all results and manipulation of the unit, during the process, by using a projector.

#### ⑤ Cables and Accessories, for normal operation.

#### ⑥ Manuals: This unit is supplied with 8 manuals: Required Services, Assembly and Installation, Interface and Control Software, Starting-up, Safety, Maintenance, Calibration & Practices Manuals.

\* References 1 to 6: ROUC + ROUC/CIB + DAB + ROUC/CCSOF + Cables and Accessories + Manuals are included in the minimum supply, enabling a normal operation.




ROUC Unit



ROUC/CIB



DAB

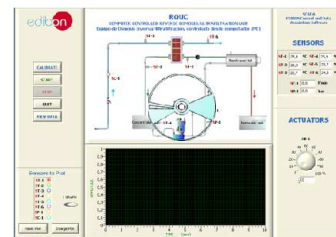

ROUC/CCSOF

**Complementary items to the standard supply**

PLC. Industrial Control using PLC (7 and 8):

**⑦ PLC-PI. PLC Module:**

Circuit diagram in the front panel.

Front panel:

**Digital inputs(X) and Digital outputs (Y) block:**

**16 Digital inputs**, activated by switches and 16 LEDs for confirmation (red).

**14 Digital outputs** (through SCSI connector) with 14 LEDs for message (green).

**Analog inputs block:**

**16 Analog inputs** (-10V. to + 10V.) (through SCSI connector).

**Analog outputs block:**

**4 Analog outputs** (-10V. to + 10V) (through SCSI connector).

**Touch screen:**

High visibility and multiple functions.

Display of a highly visible status.

Recipe function.

Bar graph function.

Flow display function.

Alarm list.

Multi language function.

True type fonts.

Back panel:

Power supply connector.

Fuse 2A.

RS-232 connector to PC.

Inside:

Power supply outputs: 24 Vdc, 12 Vdc, -12 Vdc, 12 Vdc variable.

**Panasonic PLC:**

**High-speed scan of 0.32 μsec.** for a basic instruction.

**Program capacity of 32 Ksteps**, with a sufficient comment area.

Free input AC voltage(100 to 240 VAC).

DC input: 16 (24 VDC).

Relay output: 14 (250 VA AC/2 A).

**High-speed counter.**

**Multi-point PID control.**

**Digital inputs/outputs and analog inputs/outputs Panasonic modules.**

Communication RS232 wire, to computer (PC).

**⑧ ROUC/PLC-SOF. PLC Control Software:**

For this particular unit, always included with PLC supply.



PLC-PI

**Items available on request**

**⑨ ROUC/CAL. Computer Aided Learning Software (Results Calculation and Analysis).**

**⑩ ROUC/FSS. Faults Simulation System.**



**Software Main Screens**

Main screen

The main screen displays the ROUC (Reverse Osmosis/ Ultrafiltration Unit) interface. It includes a schematic diagram of the system with components like Water Inlet, Concentrate tank, Membrane Unit, and Permeate tank. On the left, there are control buttons: CALIBRATE, START, STOP, QUIT, and VIEW DATA. Below these are 'Sensors to Plot' options for ST-1 through SC-1. A graph shows Amplitude vs. Time (secs). On the right, there are sections for SENSORS and ACTUATORS. The SENSORS section shows real-time data for ST-1 through ST-6 and SC-1, SP-1. The ACTUATORS section shows a control knob for AB-1.

**Note:** ST= Temperature sensor. SC= Flow sensor. SP= Pressure sensor. AB= Pump.

Examples of Sensors Calibration screens

The CALIBRATION window shows settings for sensor SP-1, including Gain (6,87973), Offset (0,623521), and Full Scale (50). The MULTICALIBRATE window shows a table of sensor calibration data:


Reference	Sensors	Volts	Calibrated	Err (%)
<input checked="" type="checkbox"/>	ST-1	0,2343	22,0472	1,85
<input type="checkbox"/>	ST-2	0,2366	23,9007	0
<input type="checkbox"/>	SP-1	0,0192	0,0618	23,84
<input type="checkbox"/>	SP-2	-0,0429	0,0539	23,85
<input type="checkbox"/>	SC-1	-0,0002	-0,4998	24,59
<input type="checkbox"/>	flow	0,0965	-0,0133	25,91
<input type="checkbox"/>	SC-1	-0,0698	-4,2332	28,13
<input type="checkbox"/>	SF-1	-0,0603	1,65862	22,24
<input type="checkbox"/>	ST-5	0,0013	0,008	23,89
<input type="checkbox"/>	ST-6	0,0327	3,057	20,84
<input type="checkbox"/>	ST-7	-0,0077	1,0331	22,87
<input type="checkbox"/>	SW-1	0,0417	2,876	21,03
<input type="checkbox"/>	SC-1	-0,0121	-10,6794	24,58
<input type="checkbox"/>	AN-1	0,0301	0,0301	23,87
<input type="checkbox"/>	AN-1	-0,0162	-0,0162	23,82
<input type="checkbox"/>	AN-1	0,0547	0,0547	23,85

The MULTICALIBRATE window also includes a table for GAIN, OFFSET, and p values for various sensors.


Continue...

Some typical exercises results

Pump working graph without concentrate refrigeration. The heat exchanger is switched off.

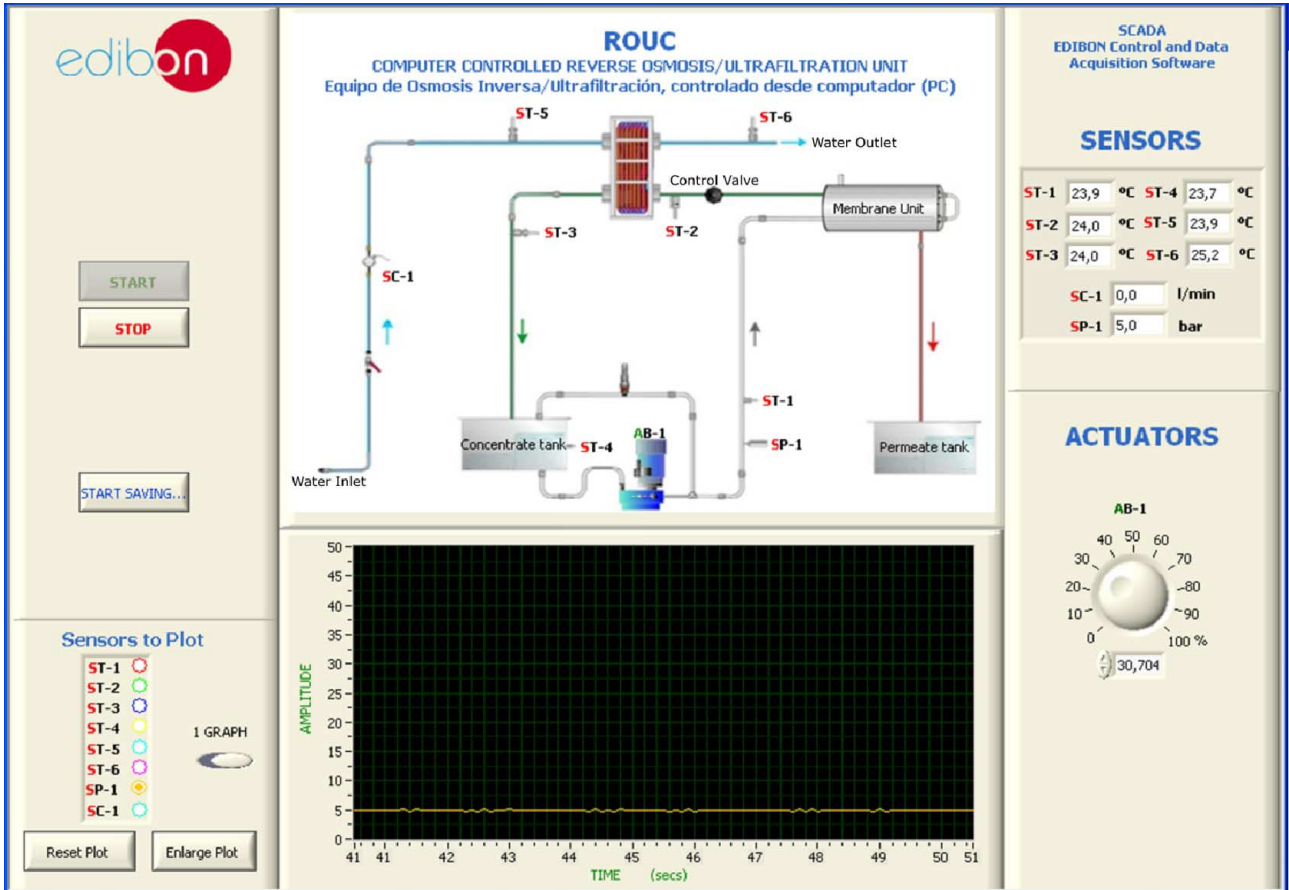

Some minutes later, after switch on the heat exchanger, the concentrate temperature is optimum for a reverse osmosis process.

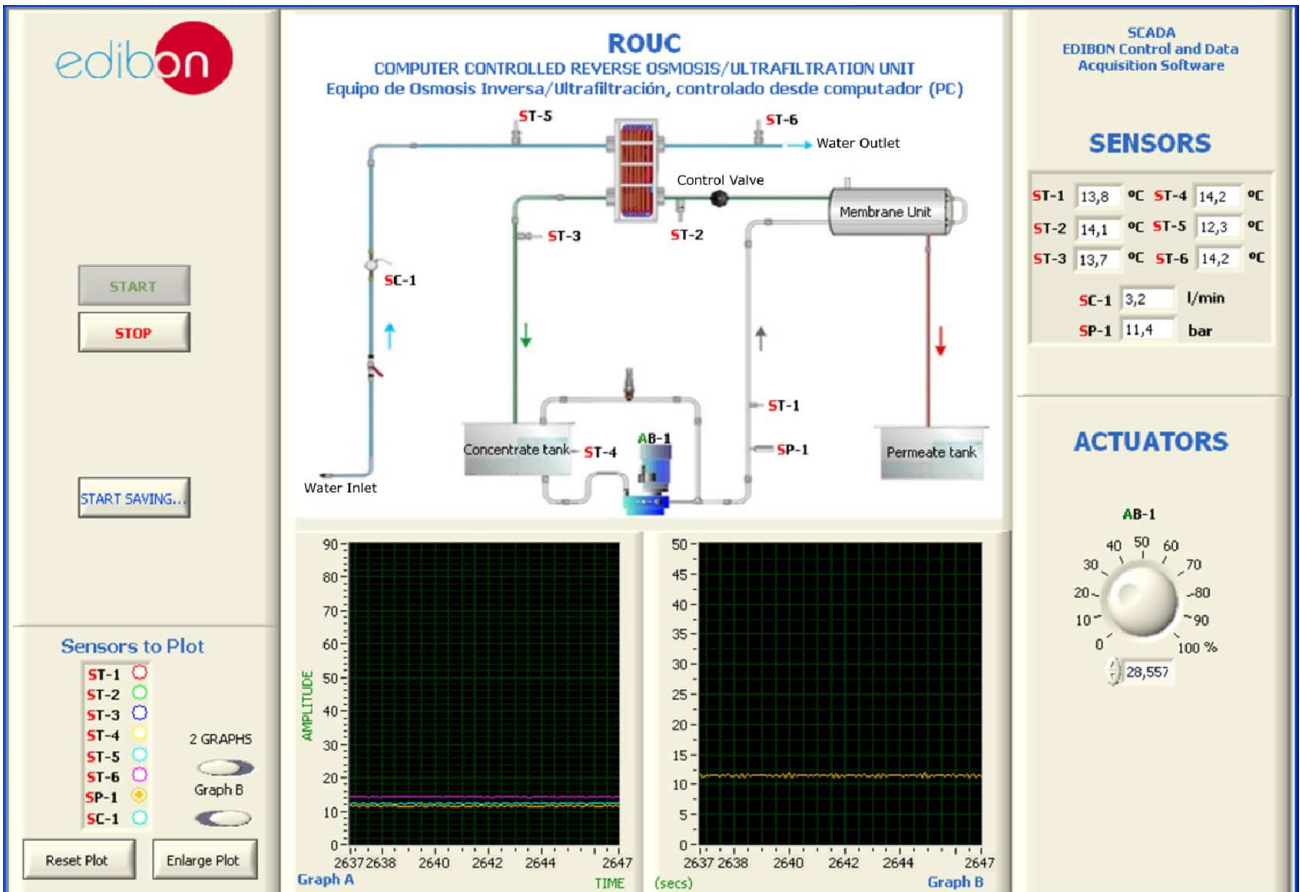


Continue...

Some typical exercises results

Membrane Module Software Protection

Besides of the protections includes in the unit, it also has software protection. If the unit reach high pressure (more than 50 bar) the software stops the process.

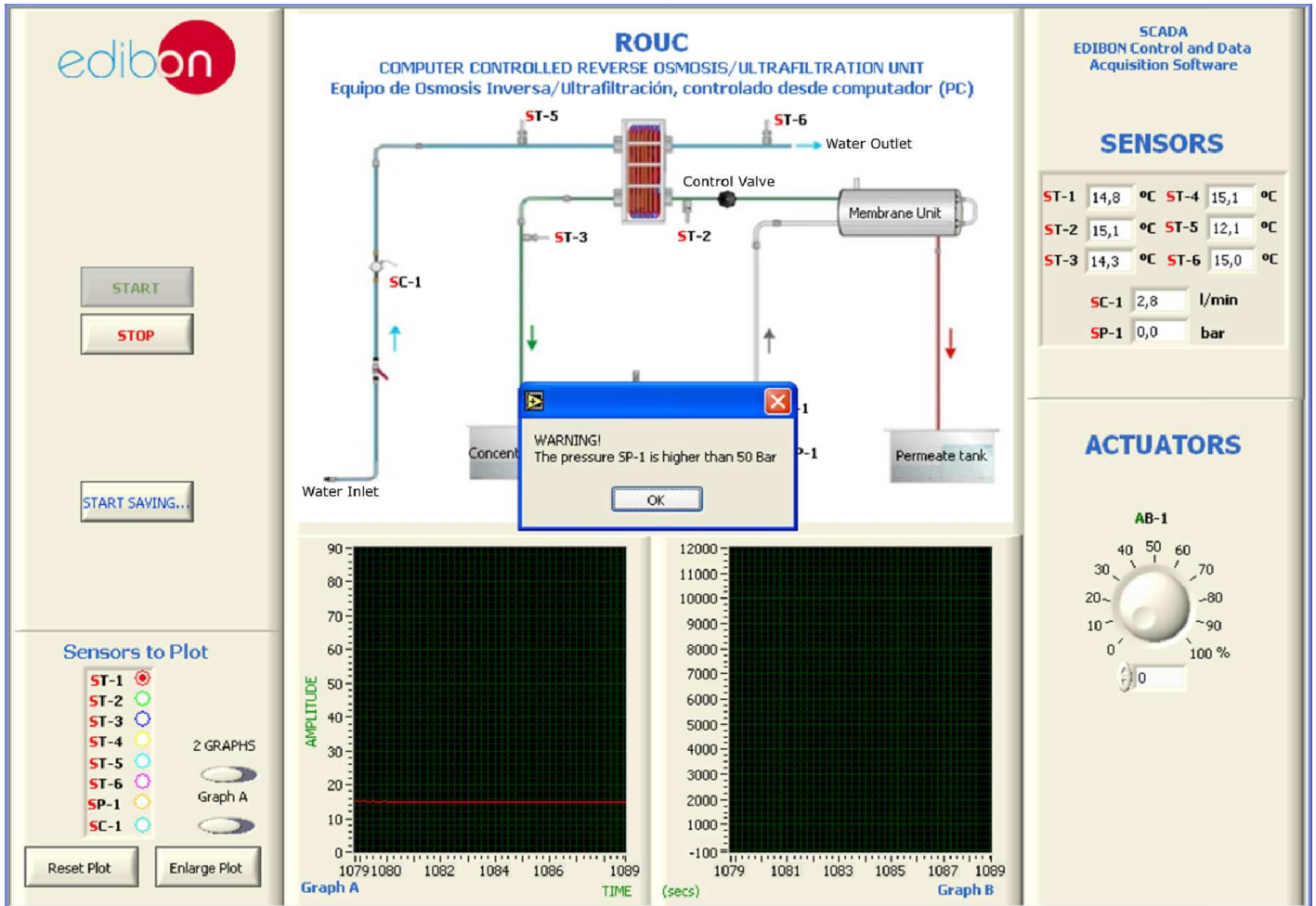







## REQUIRED SERVICES

- Electrical supply for the unit: 3-phase, 380V./50Hz or 220V./60 Hz.
- Water supply and drainage.
- Computer (PC).

## DIMENSIONS & WEIGHTS

- ROUC. Unit :
  - Dimensions: 800 x 800 x 1000 mm. approx.
  - Weight: 180 Kg. approx.
- Control Interface Box:
  - Dimensions: 490 x 330 x 310 mm. approx.
  - Weight: 10 Kg. approx.
- PLC Module (PLC-PI):
  - Dimensions: 490 x 330 x 310 mm. approx.
  - Weight: 30 Kg. approx.

\* Specifications subject to change without previous notice, due to the convenience of improvements of the product.



C/Del Agua, 14. Polígono Industrial San José de Valderas.  
28918 LĒGANĒS. (Madrid). SPAIN.  
Phone: 34-91-6199363 FAX: 34-91-6198647  
E-mail: edibon@edibon.com WEB site: www.edibon.com

Issue: ED01/10  
Date: February/2010

REPRESENTATIVE: