


## Two Pinned Arch Equipment



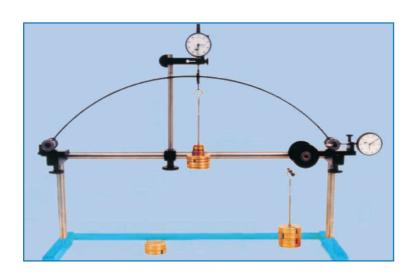

The apparatus enables the student to determine experimentally the horizontal component of the abutment thrust of a simple two hinged arch beam.

The beam is supported on ball bearing rollers attached to each end of the beam and the horizontal movement of the free end is indicated by a dial gauge so that the beam can be returned to it's original unloaded span.

The horizontal thrust force is applied to the free end of the beam by means of masses attached to a cord passing over ball bearing pulleys.

Varying loads may be applied to the beam by means of load hangers and masses and a dial gauge is provided to enable the vertical displacement to be measured.

Dimensions: 750 x 204 x 540mm

Nett WT: 9.25kg

Note: New catalogue available soon

\* Specifications subject to change without previous notice, due to continuous improvements of the product.

## **EDIBON**

**TEACHING EQUIPMENTS** 

C/San José 11-13, 28921 ALCORCON (Madrid) SPAIN Phone: 34-91-6198683 FAX: 34-91-6198647E-mail: edibon@edibon.com Web site: www.edibon.com

**DISTRIBUTED BY:** 



Edition: ED 01/01